
Panasonic

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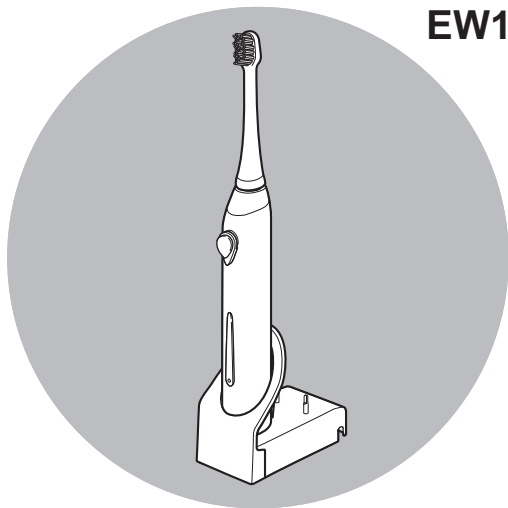Safety Instructions

Please read the safety instructions before using this product. The safety instructions provided here are to help you use this product safely and correctly, and also to prevent you from injuring yourself or others.

The safety instructions shown here are classified as follows.

WARNING: Indicates warnings that, if not heeded, may lead to death or serious injury.

CAUTION: Indicates precautions that, if not followed, may lead to serious injury or property damage.

Illustrated indications

⊘ indicates actions that you must NOT take.

● indicates actions that you should take or instructions that you should follow.

WARNING

Unplug

- Be sure to switch off the product and remove the plug from the electrical outlet before performing any maintenance. Do not touch or pull out the plug with wet hands. Failure to

do so may result in electrical shock or other injury.

Be sure to follow the instructions below.

- Clean the plug periodically so that dust does not accumulate in it. Such dust buildup may trap moisture over time and reduce the effectiveness of the insulation, thus resulting in overheating and fire.
- Use the product at the rated voltage. Failure to do so may result in accidents such as fire, electrocution and/or personal injury, including burns.

Prohibited use

- Never operate this product if it has a damaged cord or plug, if the cord or plug is warm, or if the plug fits loosely in the electrical outlet. Failure to do so may result in electrocution, short-circuit and/or fire.
- Do not damage, break or modify the cord. Do not bend, pull or twist the cord with unnecessary force. Do not place heavy objects on the cord or allow the cord to get pinched between objects. Doing so may cause fire, electrocution, short-circuit and/

or burns.

Never disassemble the product.

- Do not modify, disassemble or repair the product at any time except when disposing of it. Doing so may result in fire, electrocution and/or personal injury.

For repair or maintenance, take the product to the store where you purchased it or to an authorised service centre.

Do not use in areas where water is present.

- Do not use the charger in the bath. Doing so may result in electrocution and/or short-circuit.
- Do not wash the charger. Doing so may result in electrocution.
- Do not reach for the charger that has fallen into water. Unplug the charger immediately to avoid the risk of electrocution.

CAUTION

Be sure to follow the instructions below.

- Use this product only for its

intended use as described in this manual. Failure to do so may result in accidents and/or personal injury.

- Keep the product out of reach of infants. Failure to do so may result in burns and/or electrocution.
- Use the accessories recommended by the manufacturer. Failure to do so may result in accidents and/or personal injury.

Unplug

- Unplug when the product will not be used for more than a week. Failure to do so may result in electrocution and/or fire caused by electrical leakage as a result of inadequate insulation.

Prohibited use

- Never unplug it by pulling the cord. Doing so may result in electrocution and/or short-circuit.
- Discontinue use immediately if the product has been dropped, hit or otherwise damaged through impact. Such damage may cause electrocution and/or fire.

- This product should not be used by infants, or by those who have limited sensation inside their mouth. Doing so may result in accidents or personal injury.
- Don't let infants use. Doing so may result in personal injury.
- Do not place the toothbrush forcibly against your teeth or gums, or do not leave the toothbrush in one location in the mouth for a long time. Doing so may damage the teeth and gums. People with tooth or gum disease or people receiving oral treatment should consult their dentist before use.
- Do not brush loose dentures or teeth that are under dental treatment. Also, keep the brush away from any sores or other sensitive areas inside the mouth. Failure to observe these precautions may aggravate the condition and/or result in an accident.
- Do not use this toothbrush with toothpaste that contains a large amount of abrasives for removing stains from teeth, or with toothpaste that contains a lot of whitening agents. Failure to observe these precautions may damage your teeth or

gums.

- Do not place any metal objects, such as coins and clips, on the charger. Doing so may cause the objects to heat up, resulting in burns.
- Do not place the plastic portion of the toothbrush against the teeth or gums. Doing so may damage your teeth or gums.

For your personal dental care

Thank you for purchasing the Sonic Vibration Toothbrush.

The Sonic Vibration Toothbrush is a new kind of electric toothbrush designed for those who wish to keep their teeth clean and healthy. This product will enable you to maintain good dental hygiene, which is as important as other forms of personal hygiene.

Sonic Vibration Toothbrush Technology

The Sonic Vibration Toothbrush incorporates innovative motor technology that uses magnetic attraction and repulsion to provide outstanding performance. The brush's 26,000 strokes per minute help remove plaque gently yet

effectively.

The intensity of the toothbrush vibration is automatically weakened when the brush is pressed too hard against your teeth.

PARTS IDENTIFICATION

- Ⓐ **Toothbrush shaft**
- Ⓑ **Waterproof rubber seal**
- Ⓒ **On/off switch**
- Ⓓ **Charging indicator**
- Ⓔ **Power cord**
- Ⓕ **Charger**
- Ⓖ **V-shaped brush**
- Ⓗ **Identification ring**

BEFORE USE

This high-speed, vibrating, electric toothbrush is different from conventional electric toothbrushes.

Therefore:

- You only need to lightly press the bristles against your teeth for brushing.
- When first using this product, some may experience a tickling or tingling sensation caused by the high-speed vibrations.

The sensation should subside after you use the toothbrush several more times.

Installing the charger

Cord (Fig. 1)

Place the cord in either the right or left cord guide, depending on where the charger will be placed.

Installing near a sink or on a shelf (Fig. 2)

Place the charger on a level, secure surface. Insert the plug into an electrical outlet.

CAUTION:

- Do not use the charger in a shower stall or other area with running water that come in direct contact with the charger.
- Do not remove the plug from the electrical outlet with wet hands.
- Be sure to use only a 100-240 VAC power supply. (110-120 VAC for Taiwan)

CHARGING

Be sure to charge the main unit for at least 15 hours:

- 1) before using for the first time;
- 2) When you have not used the toothbrush for more than three months; or
- 3) if the toothbrush runs out of power during use.

Be sure to turn off the toothbrush before charging.

Place the toothbrush upright in the charger. (Fig. 3)

The charger indicator will light up red.

- The toothbrush charges whenever it is placed in the charger.
- The indicator light does not go off even after the toothbrush is fully charged.
- The toothbrush will provide about 40 minutes of use when the battery is fully charged.
- The operating duration will become shorter as the battery gets older.

CAUTION:

Do not leave metal objects, such as coins or paper clips, on the top of the charger. They may become very warm.

NOTE:

- During charging, the toothbrush will get warm, but this is normal.
- Use only the exclusive charger provided.
- If the toothbrush operates only for a short time after being fully charged, this means that the battery is nearing the end of its life cycle. Take the unit to an authorised service centre to have the battery replaced.
- Use the toothbrush at an ambient room temperature of 0° to 35°C.

HOW TO USE

1. Place the brush head firmly onto the toothbrush shaft. (Fig. 4)
 - There is a small space between the brush and the main unit. This is normal and is necessary for the head to vibrate properly.
2. Place the toothbrush in your mouth.
 - Apply small amount of normal toothpaste to the toothbrush, if desired, before turning the power on to prevent the toothpaste from spattering.

3. Turn the on/off switch to “on” position to turn on the toothbrush.
4. When you have finished brushing, turn the on/off switch to “off” position to turn off the power.
 - To prevent foam from spattering, turn the toothbrush off before removing it from your mouth.

Strong pressure preventive function

Pressing the brush bristles too hard against your teeth will activate this feature, causing the intensity of the vibration to be weakened automatically as a safeguard against excessive pressure.

PRECAUTION:

- When you first begin using the toothbrush, you may experience slight bleeding of the gums even though your gums may be healthy. This is because your gums are receiving stimulation by the toothbrush for the first time. The bleeding should stop after one or two weeks of use. Bleeding that persists for more than two weeks may be an indication of an abnormal

condition of your gums, such as pyorrhea. In this case, stop using the toothbrush and consult your dentist.

- Do not allow hair products (e. g., gel or oil) or hand lotion to come in contact with the unit as they may damage its surface.

NOTE:

Replacing the brush head every 3 months is recommended.

If the bristles are bent to one side, the toothbrush is not very effective for brushing. Replace the brush with a new one when the bristles are bent or worn.

TIPS ON USAGE **V-Shaped Brush**

Brush carefully between teeth or around teeth that are not straight. (Fig. 6)

Place the bristles of the toothbrush gently against your teeth.

Front teeth, and between teeth and gums:

- Place the bristles at 45°. (Fig. 7)

Brush the chewing surface of teeth:

- Brush at a 90° angle. (Fig. 8)

Brush the back of teeth or

where teeth are not straight:

- Brush vertically.
- Move the toothbrush to match the angle of the teeth. (Fig. 9)

AFTER USE

Toothbrush

1. Wash the toothbrush with water and shake well to remove the water. (Fig. 10)
2. Return the brush to the brush stand.
3. Use a cloth to clean the main unit. Wash hard-to-remove stains with water.
 - Do not submerge the main unit in water for a long time.
 - After washing with water, wipe dry with a cloth.
 - Do not scratch or damage the waterproof rubber seal. Doing so allow water to seep inside and damage the main unit.
4. Place the main unit on the charger.

If the identification ring comes off, put it back on. (Fig. 11)

- Do not forcibly pull the identification ring since it can break easily.
- Using the toothbrush without the identification ring has no effect on performance; the ring is simply for user identification purposes.

MAINTENANCE

CAUTION:

- Be sure to remove the plug from the electrical outlet before cleaning the product. Failure to do so may result in electrical shock or other injury.

Charger

Use a cloth to wipe off the charger.

- If cleaning requires a detergent, use a mild, neutral cleanser.
- Do not use harsh chemicals such as thinner, benzine or alcohol as these will discolour the surface.

About once every six months wipe the prongs of the power plug with a dry cloth.

CAUTION:

- Do not wash or submerge the charger in water.
- Clean the power plug periodically. Dirt and dust buildup may cause electrical shock or reduce the effectiveness of the insulation, thus resulting in overheating and a fire.
- The supply cord can not be replaced. If the cord is damaged, the appliance should be scrapped.

TROUBLESHOOTING

Problem	Possible cause	Remedy
The intensity of the vibration weakens.	The toothbrush is pressed hard against your teeth and excessive pressure preventive function is activated.	Place the brush gently against your teeth.
The toothbrush operates only for a short period of time even immediately after being charged.	The main unit was not positioned straight up in the charger and was unable to charge properly. (See page 34.)	Stand the main unit up straight so that the bottom touches the charging base, and make sure the charging indicator comes on.
	The battery has reached the end of its life cycle (about 3 years).	Take the unit to an authorised service centre to have the battery replaced. (Do not attempt to replace the battery yourself since you may damage the unit's waterproofing or internal mechanisms.)
	Charging time was insufficient.	Charge the main unit for at least 15 hours. (See page 34.)

Problem	Possible cause	Remedy
The toothbrush doesn't operate.	You have just purchased the product or it has not been used for more than 3 months.	Charge the main unit for at least 15 hours. (See page 34.)
The direction of the brush movement changes. (The brush head moves from side to side.)		Not a malfunction. Remove the brush head once and put it back to the correct position.
The tooth brush is not working properly. (For example, its movement is uneven.)		Place the main unit in the charger for at least 15 hours, then turn it on again to see if it operates properly.

If the product still does not function properly, stop using it and contact an authorised service centre.

WARNING:

Do not disassemble the product at any time other than when disposing of it. Doing so may result in short circuit, fire, or injury. For repair or maintenance (such as battery replacement), take the product to an authorised service centre.

The supply cord of the charging adaptor cannot be replaced; if the cord is damaged, the entire charger should be discarded.

DISPOSAL OF THE PRODUCT (Fig. 12)

The power source of this product is nickel-metal hydride battery which can be recycled. Please dispose of the device and the battery in accordance with national or public regulations or environmental protection or recycling regulations.

WARNING:

Do not disassemble the product at any time other than when disposing of it. Failure to observe this point may result in short circuit, fire, or injury. For repair or maintenance (such as battery replacement), take the product to an authorised service centre.

Removing the Ni-MH battery before disposing of the product

- Before removing the battery, be sure to completely drain the unit of power.
- Use a screwdriver to disassemble the unit to dispose of the battery. Remove ①, ②, ③ and ④, ⑤. Pressing the toothbrush shaft (see page 2), remove the two hooks ⑥. Remove the ⑦ and ⑧ as shown in Fig. 12.

- Be careful not to short-circuit the negative and positive poles of the battery.

SPECIFICATIONS

Rated voltage:

100-240 V AC
(110 - 120 V AC for Taiwan)
50 - 60 Hz

Battery:

1 Ni-MH battery (1.2 V DC;
1,100 mAh)

Charging time: 15 hours

Input current: 35 mA

Power consumption:

2 W (when fully charged)

Brush rotating speed:

approx. 26,000 strokes/min.
(with no load)

Brush movement width:

approx. 0.8 mm

Operating time:

approx. 40 minutes of
continuous use (with 300 g
load at a room temperature of
20°C after the battery has
been fully charged.)

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Panasonic Electric Works (Thailand) Co., Ltd.